

STIMULUS LIGHT

ENV-621L

WIRING GUIDE

DOC-221

Rev. 1.1

Copyright ©2012

All Rights Reserved

Med Associates Inc.

P.O. Box 319

St. Albans, Vermont 05478

Phone: 802.527.2343

Fax: 802.527.5095

www.med-associates.com


CHAPTER 1 | CONNECTIONS

OPERATE CABLE

Connect the Stimulus Light cable (part number SG-222) to any available Output on a Standard MED Connection Panel. When a signal is sent to the Output, the stimulus light will illuminate.

PINOUTS

Molex Pin #	Function	Wire Color
1	Not Used	Black
2	Operate	White
3	+28 Volts	Red

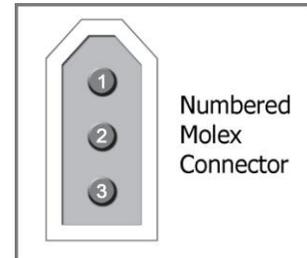

APPENDIX A | CONTACT INFORMATION

Please contact MED Associates, Inc. for information regarding any of our products.

Visit our website at www.med-associates.com for contact information.

For technical questions, email support@med-associates.com.